

MOLECULAR WEIGHT DNA MARKERS/LADDERS (100 bp)

Cat No. 110381

DATA-SHEET

DESCRIPTION

100 bp DNA Ladder is specially designed for determining the size of double strand DNA from 100 to 2,000 base pairs. The 100 bp DNA Ladder consists of 13 double strand DNA fragments ranging in sizes from 100 to 1,000 bp in 100 bp increments, and additional fragments of 1,200, 1,600, 2,000 bp. The 500, 1,000 and 2,000 bp bands are two to three times brighter for easy identification.

Supplied DNA molecular weight marker is designed to ensure easy, accurate and reliable sizing of linear double stranded DNA and posses the following features:

- Bands progress in even, easy to remember increments.
- Ladder is supplied with an optimized concentration of loading buffer.
- DNA molecular weight marker is stable for 6 months, when stored at 4°C.
- Ladder consists of ten fragments that range in size from 100 to 2,000 bp.

Storage Conditions

Store at -20°C

QUANTITY SUPPLIED; 200µl



Manufactured and Marketed By :

GENOME DIAGNOSTICS PVT. LTD.,
(An ISO 13485:2003 & ISO 9001:2000 Certified Company)

Up Mohal Naryal, Khasra No. 427, Opp. Divya Packers, Old Timber Depot Road,
Near sector 4, Parwanoo, Dist. Solan. H.P.

Email: genome24@rediffmail.com

Version : MWM04/001


www.genome-diagnostics.com

www.genomediagnostics.in

Tel: 00 91 01792 234285,, 234286

www.genomediagnostics.co.in

www.diagnosticsgenome.com

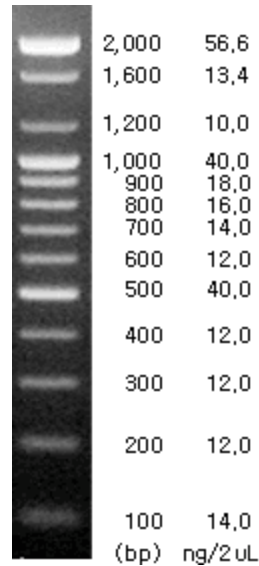

100bp DNA Ladder

Load 5–10 µl of the supplied marker into the well of Agarose or acrylamide gel.

NOTE : For research use only. Not for use in Diagnostic or therapeutic procedures.
Genome Diagnostics Pvt. Ltd., shall not in any event be liable for incidental or special damage of any kind resulting from any use except on applications.



Manufactured and Marketed By :

GENOME DIAGNOSTICS PVT. LTD.,
(An ISO 13485:2003 & ISO 9001:2000 Certified Company)

Up Mohal Naryal, Khasra No. 427, Opp. Divya Packers, Old Timber Depot Road,
Near sector 4, Parwanoo, Dist. Solan. H.P.

Email: genome24@rediffmail.com

Version : MWM04/001

www.genome-diagnostics.com

www.genomediagnostics.in

Tel.: 00 91 01792 234285,, 234286

www.genomediagnostics.co.in

www.diagnosticsgenome.com